

The *Alliana* Air Active Day Chair

Powerful mattress therapy is now mobile

Head to heel 10cm air cells alternate for deep tissue therapy

Multi-stretch PU top cover with hyperstretch coverlet for insulation and comfort



Fully adjustable day chair from total recline to upright. Side wings open out for easy transfer and cleaning (meal table optional extra - available on request)



Deep cycle gel battery for up to 2 days of therapy

Bluetooth Control Unit
Powerful yet silent with full function remote control and monitoring from any Apple or Android device

Pressure Maps Sacrum and Lower Back - 92kg Male

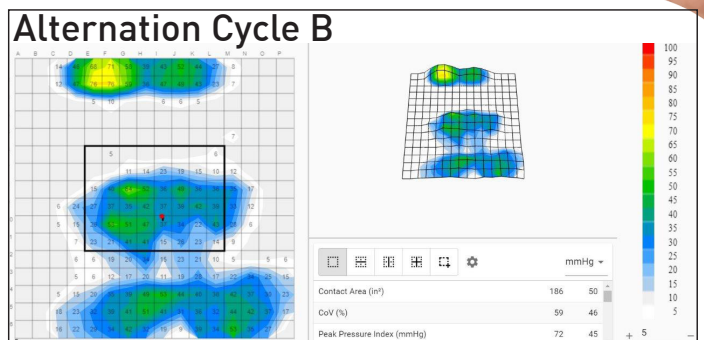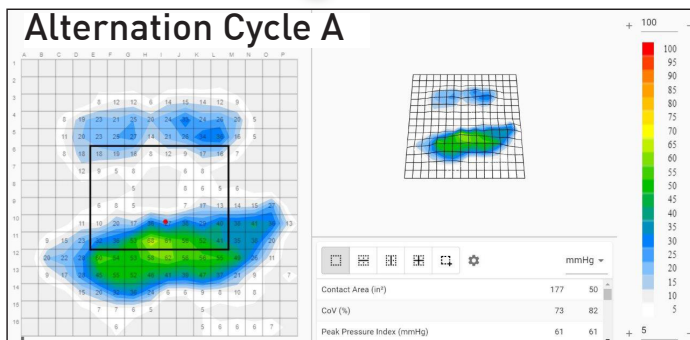

Independent tests 23/6/21

437 Nepean Highway, Brighton East VIC 3187

t: 03) 9596 6377 f: 03) 9596 6343 e: info@pegasushealthcare.com.au

www.pegasushealthcare.com.au



The Alliana Air Active Day Chair

Specifications:



Order Code: Alliana-DCS-001 Bariatric: Alliana-DCB-001

Chair and Cell System:

Safe working load: 150kg Bariatric: 200kg
Direction and brake lock castors
Overall height: 1300mm
Overall width: 710mm Bariatric: 800mm
Overall depth: 560mm
Seat level height: 700mm
Seat width: 510mm Bariatric: 600mm
Seat depth: 560mm
Seat back adjust angle: 97-136°
Seat tilt adjust angle: 0-17°
Leg rest adjust angle: 0-90°
Includes headrest and foot pads
Mealtray optional

Unit weight without battery: 55kg
Shipping dimensions: 1150x700x1300mm
Warranties: Chair and control unit: 2 years
Battery/inverter/charger: 1 year

Air alternating control system:

Features:
Therapy, static (60mins) and transfer (20 mins) modes with safety timeouts to therapy.
Selectable alternating times from 10 minutes (default) to 15, 20 and 25 minutes.
Control panel auto lock mode to restrict unauthorised users.
Audio and visual alarms for leak, power and tilt sensor.
Easy bluetooth connection for remote control and monitoring with a free downloadable App.
Bluetooth charging control and monitoring with free downloadable App.
Up to 2 days charge life dependant on use and age of battery.
Easy wall charging system when not in use.
Heavy duty rear equipment bag for protection and storage



437 Nepean Highway, Brighton East VIC 3187

t: 03) 9596 6377 f: 03) 9596 6343 e: info@pegasushealthcare.com.au

www.pegasushealthcare.com.au

